


KBr - Table Press

Compact designed KBr-Table Press for pellet making is used for transmission analysis which is

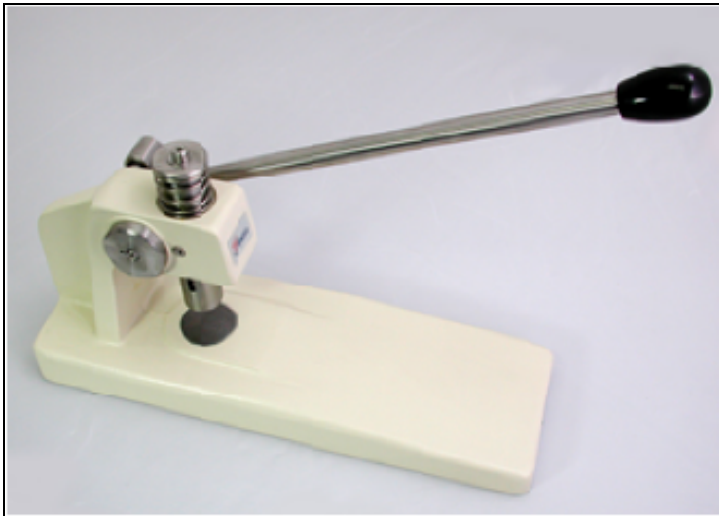

Sturdy, compact and easy To Use. The smooth action of the handle permits pellets to be pressed effortlessly with one hand and minimal pressure. Several different Die Sets (1mm, 3mm, 7mm) are available for use with the Table Press. The pellets are pressed using light operator pressure in a collar/die that can be transferred onto a pellet holder for transmission analysis.


The components that contact the die set are stainless steel Uniform and highly reproducible pellets Easy to clean and to maintain.

The QUICK HOLDER is especially designed to hold pellets in place. It uses a magnetic plate to secure the dies at the right position.

The Quick Holder was specifically designed to hold collars for 1mm, 3mm, 7mm, and 13mm pellets. It has a 2"x3" backplate, which fits into the sample compartment of all standard FTIR spectrometers.

- The magnetic back plate holds the sample in place. Easy positioning of the collar on a removable mount with a half moon cut-out.
- Includes 2 interchangeable mounts of different width.
- An opening at the upper end of the holder is used for 13 mm pellets.



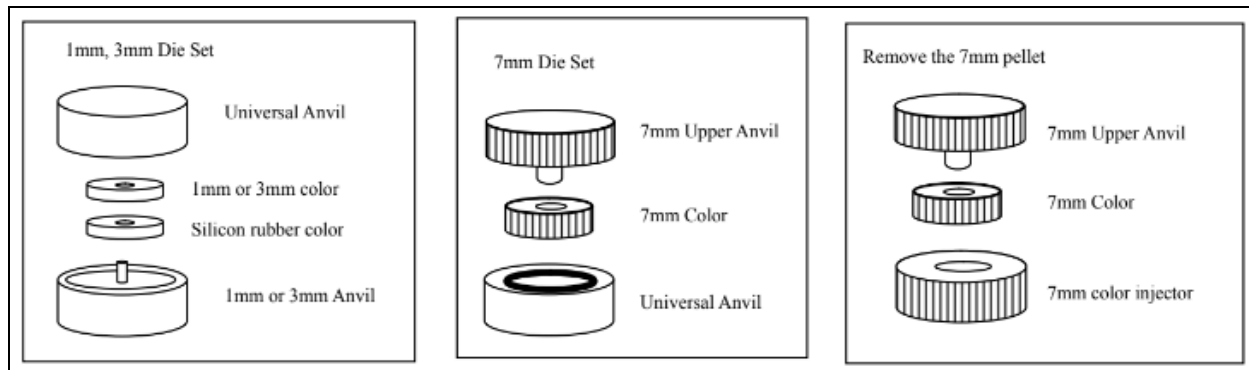

ORDERING INFORMATION

Part Number

STJ-0130
 STJ-0161
 STJ-0163
 STJ-0164
 STJ-0165
 STJ-0156

Description

Table Press Alone
 3 mm Die Set
 7 mm Die Set
 3-7 mm Die Set
 3mm Silicon Collars (6 pieces/package)
 Quick Holder