

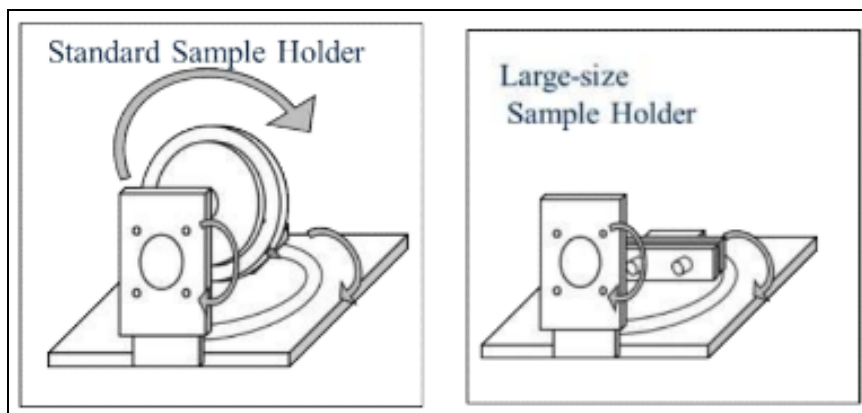
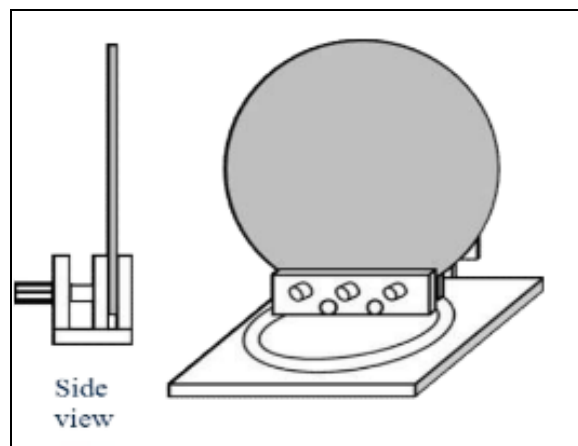
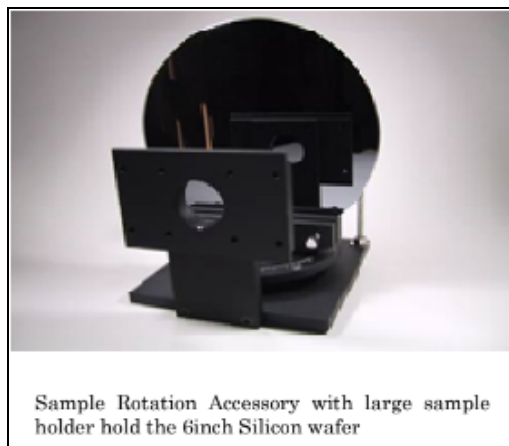


Sample Rotational Accessory

Samples can be rotated in two directions. This is a powerful tool for making transmission measurements on samples with interference fringes. This accessory can be used with polarizer for molecular orientation studies. This accessory has large size sample holder which can be hold the large sample such as silicon wafer.

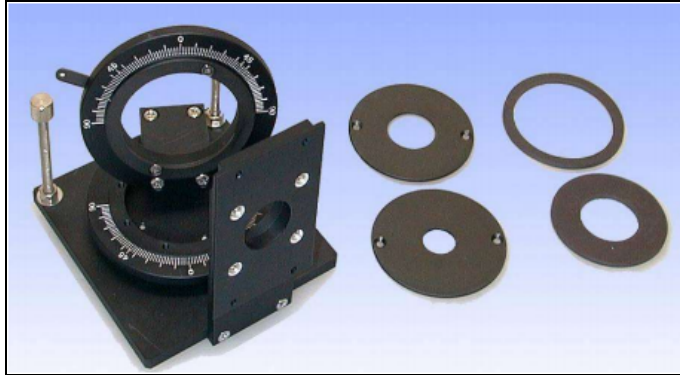

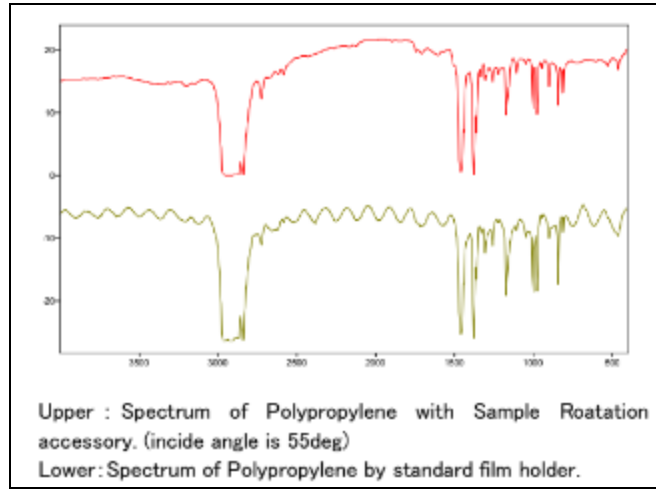

Specifications:

Back plate size : ϕ 50,24,15mm and 10
Magnet holder size : ϕ 52,28 and 18mm (T-2mm)
Scale : Each 2degree

ORDERING INFORMATION

Part Number	Description
STJ-0141	Sample Rotation Accessory