

High-Performance Grinder Mixer

Super ShakIR works well for making the pellet of the sample with KBr chunk and powder.




SuperShakIR is a Grinder mixer which produce fine powder. It is able to crush the sample, KBr Chunk, powder and mixing them for making the transparency pellet. It can also make fine powder from solid sample because of capability of high power grind. Fine powder is good for measurement in ATR method and it is possible in this device.

Features:

- Easy-Handling & Reasonable price.
- One micro tube construction can crush sample perfectly even KBr chunk.
- Sounds from Super ShakIR is very low compared to others.



Grinder Mixture Super ShakIR

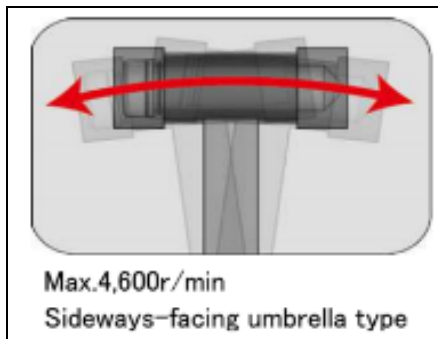

Specifications:

Speed setting:	2,500, 3,000, 4,000, 4,200, a 4,600r/min
Concussion method:	Sideways-facing umbrella type concussion method
Timer setting:	Six setting (5, 10, 15, 30, 45, 60 seconds)
Sample Tube:	NEJIRO Micro Tube: 1.5/2.0ml
Bead:	SUS

ORDERING INFORMATION

Part Number	Description
SS-01-S	SuperShekIR Including Super ShakIR, SUS tube EP, SUS Bowl bead (50pcs), SUS bullet-shaped bead EP