



## Micro Vice SliceIR – Easy Cross Sectioning:


The SliceIR is a unique attachment for the MicroVice sample holder for enhanced cross section cutting. With additional available adaptors, you can easily hold your samples and enlarge the

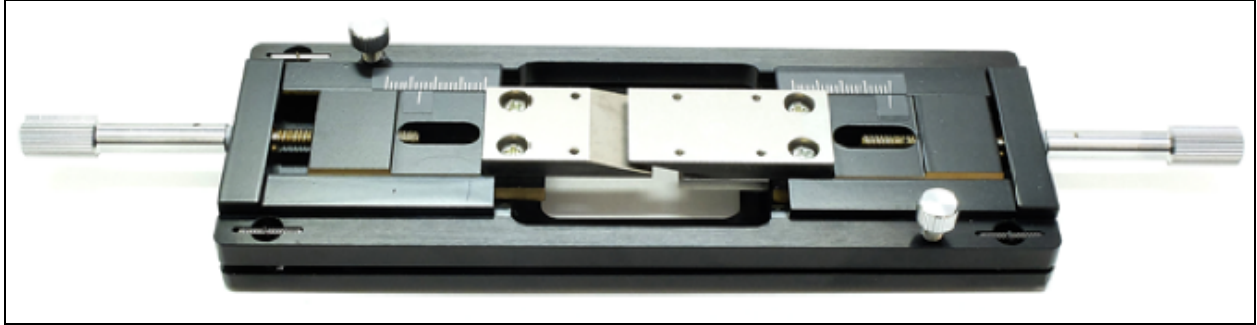

cutting surface of your sample. The other benefits are:

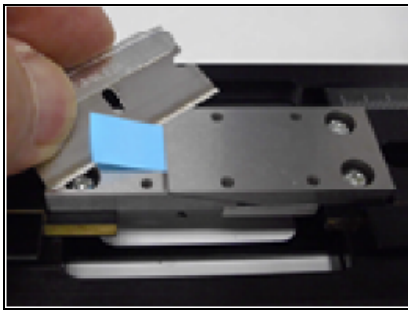

- Holding a sample with the standard MicroVice clamps provides a cutting angle of 90°.
- Holding the sample with the SliceIR clamps changes the cutting angle to 15° or 30°.
- With the 15°/30° SliceIR clamps you will have a two/four times larger cutting surface than with the standard 90° clamps.

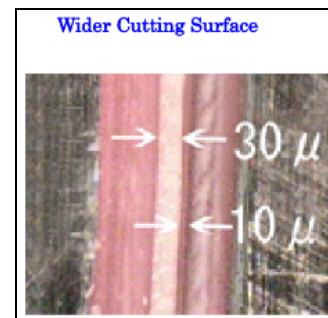

► See how it works on our  channel: <https://www.youtube.com/user/welovespectroscopy>

## ORDERING INFORMATION

Part Number	Description
STJ-0116-SL	MicroVice SliceIR (with both 30 & 15° adapter, stopper)
STJ-0116-SL30	SliceIR with 30° adapter only
STJ-0116-SL15	SliceIR with 15° adapter only
STJ-0116N	Teflon Knife (100 pieces)

